

SCOOPEXL™

Always with the patient, from the immobilization, during the transport to the emergency room up to the inside of the diagnostic and emergency departments.

USER-FRIENDLY

It allows the operator to carry out all the handling manoeuvres required for the traumatised patient in any environment. It allows to immobilize a patient on any type of ground through the intervention of only two operators. Using SCOOPEXL™ is easy, fast and safe.

COMFORT


The concave shape is perfect for the patient comfort and guarantees more contact points between the body and other immobilisation aid.

SAFETY 1

The patented Twin Safety Lock (TSL) system prevents accidental opening of the SCOOP EXL™ stretcher. SCOOPEXL™ is easy to sanitise as all its components are easy to reach.

ERGONOMICS 3

Designed and manufactured according to the most advanced ergonomics criteria, it guarantees the operator to maintain a correct posture during all rescue and transport manoeuvres of the patient.

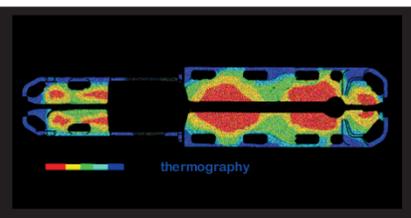

THERMAL INSULATION

SCOOPEXL™ guarantees immediate and effective thermal insulation thanks to its special plastic material of which it is made.



DIAGNOSTICS

SCOOPEXL™ can be used during X-ray imaging examinations necessary for the urgent treatment of traumatized patients. The absence of ferromagnetic materials allows its use in Magnetic Resonance Imaging.



Compatible with the 365-E immobiliser



CROSSOVER TRIALS

Studies carried out at the Center for Emergency Medicine at the University of Pittsburgh show that the degree of spinal immobilization with SCOOPEXL™ stretcher is equivalent to that obtained with a spinal board. They also show that a reduced number of manoeuvres by the operator is sufficient, thus causing less stress to the patient's backbone.

AVAILABLE COLOURS

YELLOW, RED, GREEN AND GREY



EASY TO SANITISE

COMPACTNESS AND TRANSPORTABILITY 2

Stretcher adaptability to any patient thanks to the possibility of length adjustment up to a maximum of 2010 mm. Easy to close and transport thanks to its reduced weight and minimum overall dimensions. It can be stored in an ambulance with minimum space requirements.



COMPLIANCE

Compliant with Locking systems compliant with

Not magnetic.

UNI EN 1865
UNI EN 1789



The SCOOPEXL™ stretcher is classified as a SAFETY 4 YOUR BACK device. The SAFETY 4 YOUR BACK devices are specifically designed to protect the healthcare operator's backbone.



MAIN FEATURES

- ▶ Made of special aluminium alloy and plastic polymers, HDPE.
- ▶ Patented Twin Safety Lock (TSL) system.
- ▶ Ergonomic.
- ▶ Foldable and divisible into 2 parts.
- ▶ Adjustable in length up to 2010 mm.
- ▶ Compatible with X-ray diagnostic tests and MRI.
- ▶ Easy to sanitise.
- ▶ Compatible with quick-release restraints.
- ▶ Available colours: yellow, red, green and grey.

TECHNICAL SPECIFICATIONS

Maximum length	2010 mm / 79.13 in
Minimum length	1647 mm / 64.82 in
Maximum width	432 mm / 17 in
Dimensions when folded	1200x430x80mm / 47x17x3.15 in
Weight	8 kg / 17.63 lb
Load capacity	227 kg / 500.44 lbs

EMERGENCY KIT

- 1 Bag/backpack
- 1 SCOPEXL™ stretcher
- 1 Head immobiliser 365-E
- 1 WizLoc neck brace
- 4 Snap-hook restraints 436i
- 1 Disposable thermal blanket

part no. 65EXL-EK



COMPLETE SCOOP EXL KIT

- 1 SCOPEXL™ stretcher
- 1 Head immobiliser 365-E
- 1 Bag for 365-E
- 4 Snap-hook restraints 436i

part no. KIT-65EXL/773



ACCESSORIES



Head immobiliser 365-E

Dedicated head immobiliser that guarantees the use of the stretcher as an efficient immobilization system.

part no. 365-E



Restraints

For patient safety. Thanks to the quick-release system, they reduce the application time.

part no. 0314122



Head immobiliser support

Designed to be used with the TSL expander system, it provides excellent head support during transport.

part no. 10-9900-008



Locking system

Made and certified to fix the SCOPEXL™ stretcher in horizontal or vertical position. Certified according to UNI EN 1789.

part no. 561-V / part no. 561-O

SCOPEXL™



FERNO TSL EXPANDER KIT

A high-strength aluminium alloy accessory to be applied to the Scoop EXL atraumatic stretcher. Designed to handle and transport a corpulent patient. With the Easy Lift rescue bridles, the Scoop EXL stretcher and the TSL Expander Kit can also be used with lifting devices.

MAIN FEATURES

- load capacity of 227 Kg
- assembly in a few seconds
- patented closing system to prevent accidental opening of the stretcher
- totally decontaminable
- self-releasing high-load metal restraints
- easy to store



THE KIT CONTAINS:

- 2 red expanders
- 3 restraint extensions
- 1 containment bag



for Italy

P +39 051 6860028

✉ info@ferno.it

www.ferno.it

Numero Verde
800 501711



for Schweiz

P +41-41-2596000

✉ info@ferno-schweiz.ch

www.ferno-schweiz.ch