



• Optional accessories shown

Features

- FAST Chair (ordering number 60-4120-001) can be used as a stair chair, evacuation chair or transport chair
- Ideal and cost effective solution with Ferno EZ Glide deployed track system
- Stair track can be easily attached to or removed from the chair (to minimize the weight that needs to be carried by the operator)
- All-aluminium frame construction can carry patients of up to 250 kg weight
- Two 125 mm front swivel wheels, and two 150 mm rear fixed wheels with locks
- Wheels have minimum debris collection while maintaining good traction and comfort, and provide very good maneuverability in confined spaces such as stair landings or tight corridors
- Head support and three restraints included
- Rear top handle is extendable and adjustable in four positions (1000 to 1580 mm)
- Front telescoping handles adjustable in five positions (maximum extension of 340 mm)
- Vinyl cover is easy to clean or replace
- Optional Comfort Footbar, 4-point Restraint, Mesh Basket, Rear Support Handles, and Secure Lock Mount available
- FAST chair can be also ordered without stair-track (ordering number for a chair without stair-track is 60-4100-001), and the stair-track (manual or powered) can be purchased separately later

Specifications

Height min	Overall Width	Overall Length	Seat Width	Folded L,W, Thickness	Weight without track	Weight with track	Load Limit
1000 mm	516 mm	695 mm	415 mm	1000x516x215 mm	11,2 kg	15,0 kg	250 kg

Accessories

FAST Chair Wall Mount - designed to securely hold the FAST Chair in both a static or dynamic application. Single hand release allows for quick chair fastening or detachment.

60-4190-001



4- Point Restraint for enhanced safety.

62-4400-001



Rear Support Handles provide additional gripping points for the operator when lifting the chair or maneuvering over obstacles.

65-4122-001



Comfort Footbar increases patient's comfort and sense of security.

65-4130-001



Head Support – Included Free of Charge

65-4126-001




Foldable armrests for enhanced comfort.

65-4128-001



Foot Restraint– Included Free of Charge

65-0209-002

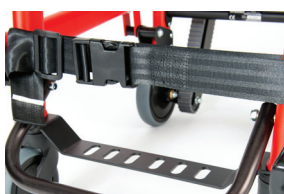

Mesh Basket for storing accessories such as oxygen bottle.

65-4123-001



Manual stair-track for a chair purchased without a stair-track.

65-0402-001



Powertraxx FAST - intelligent stair-track battery powered system that glides up to 225 kg loads up or down stairs without carrying or lifting, with up to 20 flights of stairs per charge. Chair is easy to operate when climbing or descending stairs and intuitive, easy-to-use paddle-style controls keep you in control.

65-0403-001



FAST chair has been designed to comply with and exceed international standards such as EN 1865-4 2012.

